

## Surgeon Mask Tie 05

### CONSTRUCTION

UPPER PLY : NON WOVEN + PP SPUNBOUND 18GSM

MIDDLE PLY : NON WOVEN + PP SPUNBOUND + MELTBLOWN 18GSM

LOWER PLY : NON WOVEN + PP SPUNBOUND 25GSM

LACE : NON WOVEN + PP SPUNBOUND 32GSM

### MEASUREMENT

LENGTH: 18CM ± 5%

WIDTH : 9CM ± 5% (NOT LESS THAN 16CM WHEN FULLY EXTENDED)

NUMBER OF PLEATS : 3

LACE 91CM ± 5%

NOSE STRIP : NOT LESS THAN 100MM X 2MM

## Surgeon Mask Loop 05

### CONSTRUCTION

UPPER PLY: NON WOVEN + PP SPUNBOUND 18GSM

MIDDLE PLY : NON WOVEN + PP SPUNBOUND + MELTBLOWN 18GSM

LOWER PLY : NON WOVEN + PP SPUNBOUND 25GSM

LACE : NON WOVEN + PP SPUNBOUND 32GSM

### MEASUREMENT

LENGTH : 18CM ± 5%

WIDTH : 9CM ± 5% (NOT LESS THAN 16CM WHEN FULLY EXTENDED)

NUMBER OF PLEATS : 3

LOOP : 14.5CM ±5%

NOSE STRIP : NOT LESS THAN 100MM x 2MM

## Surgeon Mask Tie 04

### CONSTRUCTION

UPPER PLY : NON WOVEN + PP SPUNBOUND 18GSM

MIDDLE PLY : NON WOVEN + PP SPUNBOUND 06GSM / 18GSM

LOWER PLY : NON WOVEN + PP SPUNBOUND 25GSM

LACE : NON WOVEN + PP SPUNBOUND 32GSM

### MEASUREMENT

LENGTH : 18CM ± 5%

WIDTH : 9CM ± 5% (NOT LESS THAN 16CM WHEN FULLY EXTENDED)

NUMBER OF PLEATS : 3

LACE : 91CM ±5%

NOSE STRIP : NOT LESS THAN 100MM x 2MM

## Surgeon Mask Loop 04

### CONSTRUCTION

UPPER PLY : NON WOVEN + PP SPUNBOUND 18GSM

MIDDLE PLY : NON WOVEN + PP SPUNBOUND 06GSM / 18GSM

LOWER PLY : NON WOVEN + PP SPUNBOUND 25GSM

LOOP : NON WOVEN + PP SPUNBOUND

### MEASUREMENT

LENGTH : 18CM ± 5%

WIDTH : 9CM ± 5% (NOT LESS THAN 16CM WHEN FULLY EXTENDED)

NUMBER OF PLEATS : 3

LOOP: 14.5CM ±5%

NOSE STRIP : NOT LESS THAN 100MM x 2MM



Surgeon Mask Loop



Surgeon Mask Tie

## HOW TO WEAR A MASK



Wear it with the coloured side facing out and white side facing in.



Ensure it covers the nose and mouth fully. Tie/loop the straps around your head or over the ear.



Use two fingers to press the concealed metal wire down to fit the shape of the nose.



If the mask is properly used, avoid touching the mask. If you want to touch the mask, make sure to wash your hands first.

## NON WOVEN SURGICAL FACE MASKS AVAILABLE IN THREE COLOURS.




Mainly manufacturing in Blue, White and Green. For customisations in colour please contact us.

## Features

We have managed to design and ensure a significantly qualitative mask for your safety. Here are the main features of our surgical masks.

- All our masks are made using an ultrasonic weld to prevent fabric de-lamination.
- All our masks have a nose strip.
- Composed of latex free, fibreglass free and hypoallergenic material.
- Standardised ear-loop for optimum periphery seal. Ample in length yet soft and strong for maximum comfort.
- Inner facing layer non-Woven fabric with Fluid absorbent property. Void of dyes, chemical and inks. Soft and odourless.
- Outer facing layer non-woven fabric with fluid repellent property. Odourless and free of lint.
- Nose Strip is malleable flat aluminium or polyethylene covered with steal wire further covered with non - woven fabric for longer wearability and to prevent it from protruding out.
- We are ISO and CE certified.

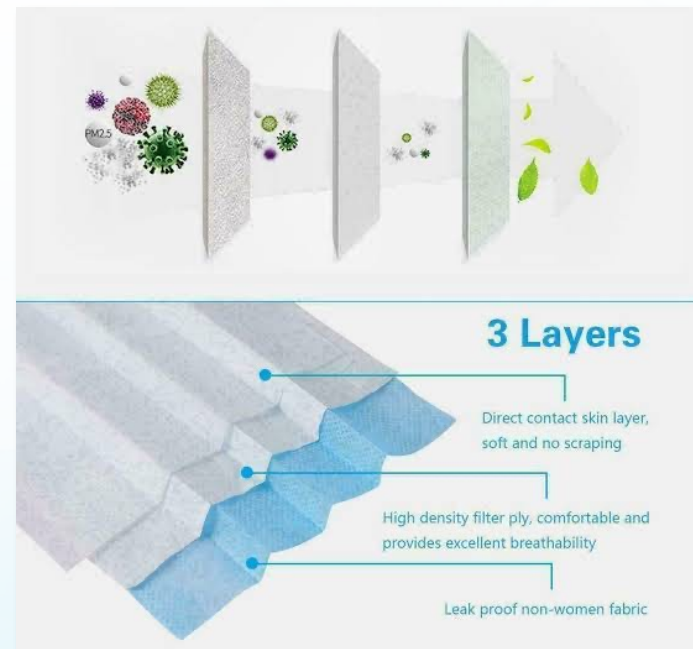

## About Us

S4 Healthcare Pvt. Ltd. an ISO 9001:2000 certified, manufacturer of disposable healthcare products. It has been a leading manufacturer, exporter and supplier since 2006 for an exclusive collection of products. Owing to the compatible and efficient use of our offerings, we have developed a trusted image in the market under the leadership of our highly experienced and skilled team. S4 Healthcare Products are designed to provide the maximum protection and comfort for its clients especially healthcare professionals and patients. We provide for flexibility in design and size based on the input of the customer. They are manufactured with utmost care under stringent conditions, adhering to ISO 9001:2000 ensuring total quality products. With a great client network of Hospitals, Corporate Customers & dealership network the company is poised to introduce many more products in the healthcare segment.

S4 Healthcare Pvt. Ltd.

23-25, Sheetal Supreme Industrial Estate No. 2 Behind Dhuri Industries No. 2 Sativali Road, Vasai (East) Thane District -401208

Tel : +91-22-40167000

Cell : +91-9323623557 and +91-9870058585

Fax : +91-22-40167001

Email : [info@S4healthcare.com](mailto:info@S4healthcare.com)

Website: [www.s4healthcare.com](http://www.s4healthcare.com)